AI-07.01 Co	ntrolserver network settings for eBox	CARSTENS Page 1 of 3
AG501	Controlserver	

Configure controlserver network settings for eBox-connection

Step-by-step tutorial

1. Connect the controlserver to the AG501 via the yellow LAN-cable (Figure 1).



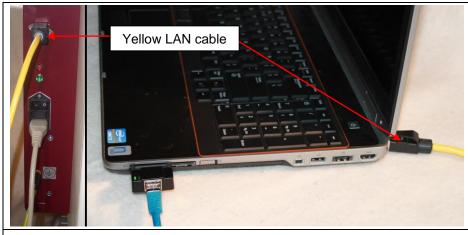


Figure 1: E-Box and Control Server cable connections

2. Right-click on the LAN-icon in the top right corner (Figure 2-1) and select *Edit Connections* in the drop-down menu (Figure 2-2).

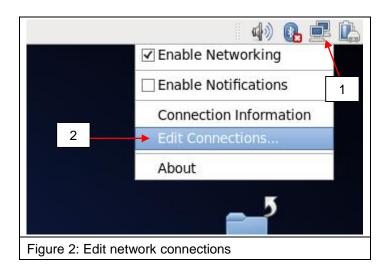


3. Select System em1 (Figure 3-1) and click on Edit (Figure 3-2).

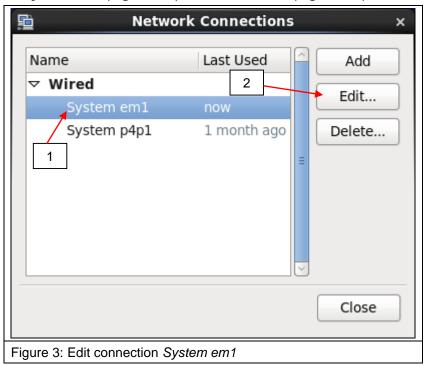


4. Select the *IPv4 Settings* tab (Figure 4-1) and configure the settings as shown in Figure 4.

Click *Routes* (Figure 4-2) to configure the routing.

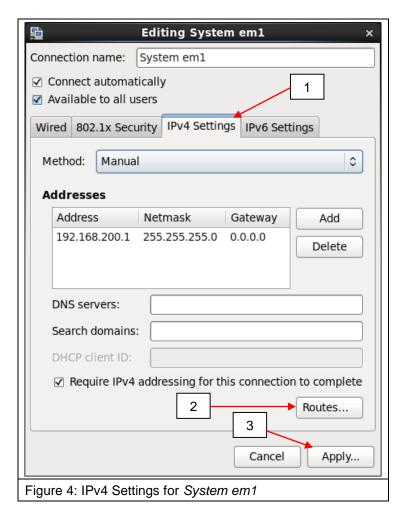

Page 2 of 3 September 21st, 2015

5. Check the *Use this connection...* as shown in *Figure 5* and click *OK*.

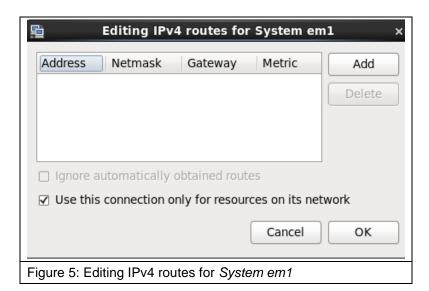

6. Click *Apply* (Figure 4-3) and type in the controlserver root-password to authenticate (Figure 6). The new settings will be applied.



7. Reboot the controlserver.

Page 3 of 3 September 21st, 2015